

## Bravo Mini

24-Hour Ambulatory  
Blood Pressure Monitoring System



### SunTech Bravo Mini

The Bravo Mini system from SunTech Medical® sets the standard in 24-hour ambulatory blood pressure monitoring (ABPM).

### Clinically Validated

The Bravo Mini system is independently tested to meet the accuracy and performance requirements of the British Hypertension Society, the European Society of Hypertension International and the ISO 81060-2 standards.

### SunTech Motion Tolerance


The Bravo Mini utilizes SunTech proprietary motion-tolerant technology, obtaining reliable readings even during motion. The Bravo Mini is validated to the ISO 81060-2 ergometry protocol to ensure accuracy and dependability.

### Smart Inflate

Smart inflate algorithm that bases measurement inflation on previous measurements to protect against over-inflation and provide superior patient comfort.

### Orbit™ Blood Pressure Cuffs

SunTech's patented Orbit cuff is specifically engineered for ambulatory blood pressure monitoring applications. This cuff is designed with a soft sleeve that gently hugs the arm to maintain cuff placement and provide the ultimate in patient comfort for an extended period of time.

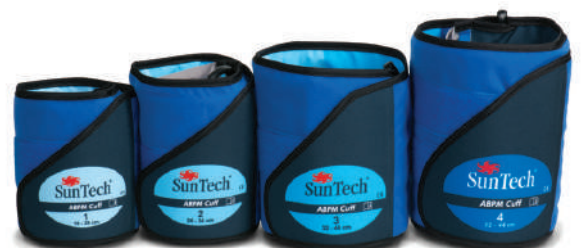

## AccuWin™ 4 SE Software

AccuWin 4 SE is a user-friendly, Windows® based program that allows for configuration, analysis, interpretation and reporting of ABPM studies.



## Visual Blood Pressure Results

The unique AccuWin™ report makes communication with patients simple and effective with easy red-green graphical representation of blood pressure study data.

## Automatic Interpretive Summary

Rapidly decipher data, creating interpretive summaries, according to AHA, ESH and Pediatric guidelines, for:

- ✓ Automated analysis of ABP levels for normality.
- ✓ White coat syndrome.
- ✓ Overnight dipper status.

## Client - Server system

Access patient information from across your network with patient data securely stored in a networked server.

## Technical Specifications

### The Bravo Mini System Includes

- ✓ Bravo Mini monitor.
- ✓ USB drive with AccuWin™ 4 SE (AW4 SE) software, quick start guide and user manual.
- ✓ 2 Orbit™ ABPM cuffs (Size 2 and Size 3) carrying pouch & patient belt.
- ✓ USB interface cable.
- ✓ Start Here Insert.
- ✓ Standard 2 year warranty.

## General Specifications

Validations:	BHS (A/A), ESH, ISO 81060-2:2013
Accuracy:	Heart Rate accurate within +/-2% or +/-3 bpm, whichever is greater. Transducer Static Pressure: ±3 mmHg
BP Technique:	Oscillometry with step deflation
BP Range:	25-260 mmHg
HR Range:	40-200 bpm
Sampling Periods:	4 adjustable (awake, asleep, and 2 special)
Time Intervals:	5, 10, 15, 20, 30, 45, 60, 90 and 120 min.
Data Storage:	Flash memory stores up to 250 readings
Data Transfer:	USB 3.0 (USB-C)
Power Source:	Two AA batteries
Weight:	233 grams, including batteries
Dimensions:	100mm x 70mm x 30 mm
PC Requirements:	USB port, Windows® 7 or later
Operating Conditions:	10°C (50) to 50°C (122°F) 20-95% RH non-condensing
Shipping/Storage Conditions:	-20°C (-4°F) to +70°C (+158°F) 15-95% RH non-condensing